




## KeriFlex<sup>®</sup>

Silicone Finger Joint Implants


PIP - MCP



Fully committed to **Hand & Wrist** only  
[www.kerimedical.com](http://www.kerimedical.com)

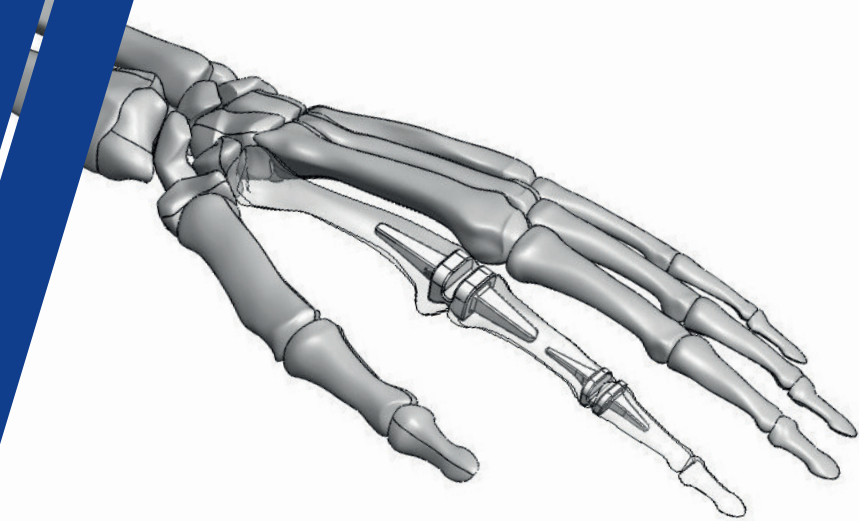

Test lab

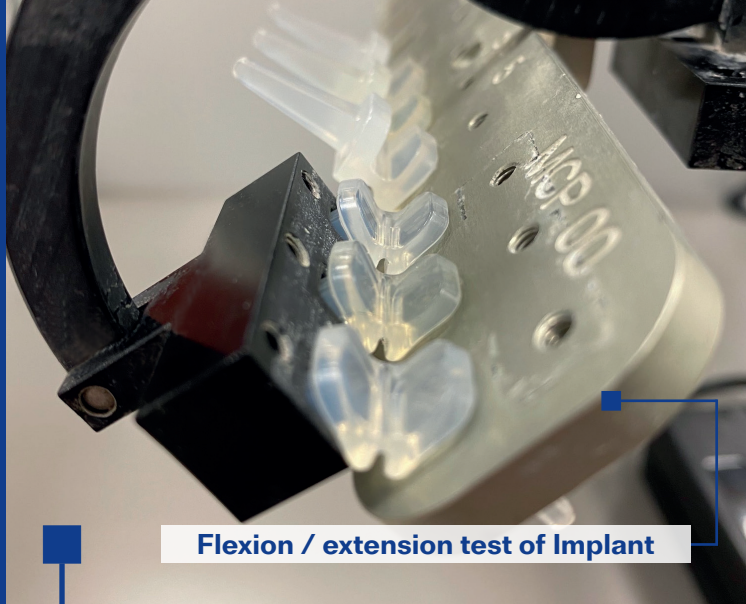

Flexion / extension test of Implant


DESIGN & RESISTANCE VERIFIED

- Full flexion / extension (0° - 90°)
- Verified resistance up to 10 million cycles (test lab)
- Respect of joint angulation at rest :  
PIP 15° and MCP 30°

- + COMFORT
- + STABILITY

KeriFlex® implants provide strength,  
flexibility and durability

MEDICAL SILICONE  
and BIOCOMPATIBLE

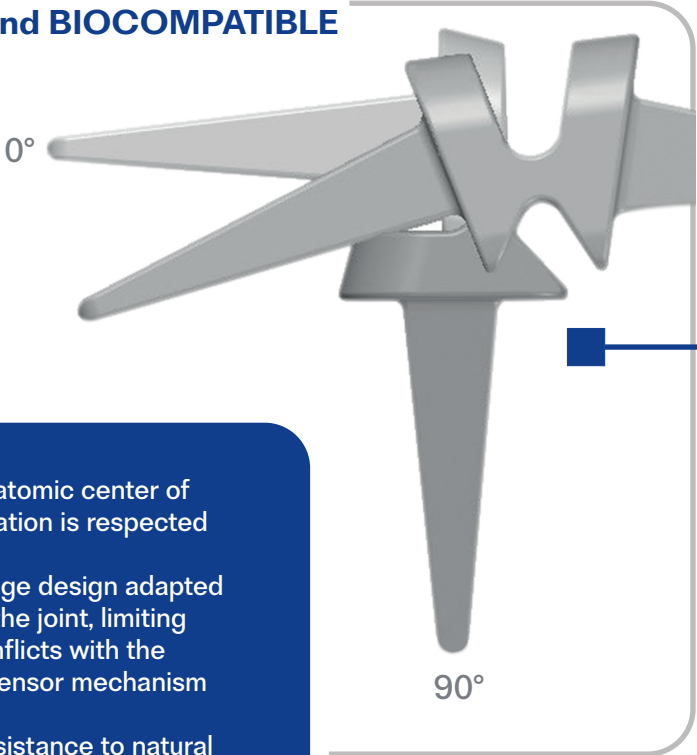

- Anatomic center of rotation is respected
- Hinge design adapted to the joint, limiting conflicts with the extensor mechanism
- Resistance to natural compression stress of the joint

The new generation  
of Finger Joint Implants

Medio-lateral stability  
and joint mobility


High performance  
medical silicone  
for high mechanical  
resistance  
to physiological efforts

Complete size range to  
meet all anatomical needs

12 SIZES AVAILABLE


PIP

5 sizes of implants

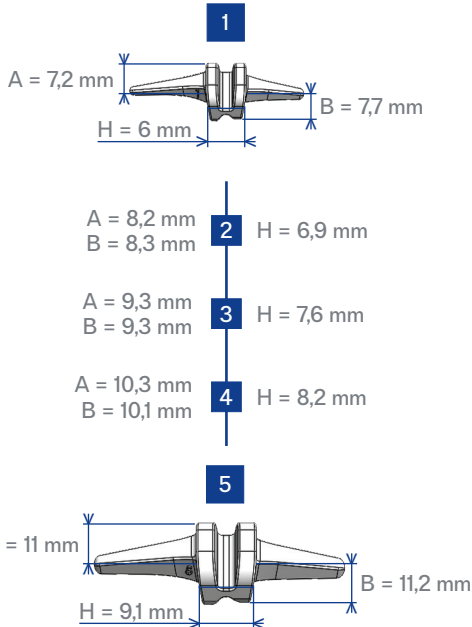

MCP

7 sizes of implants

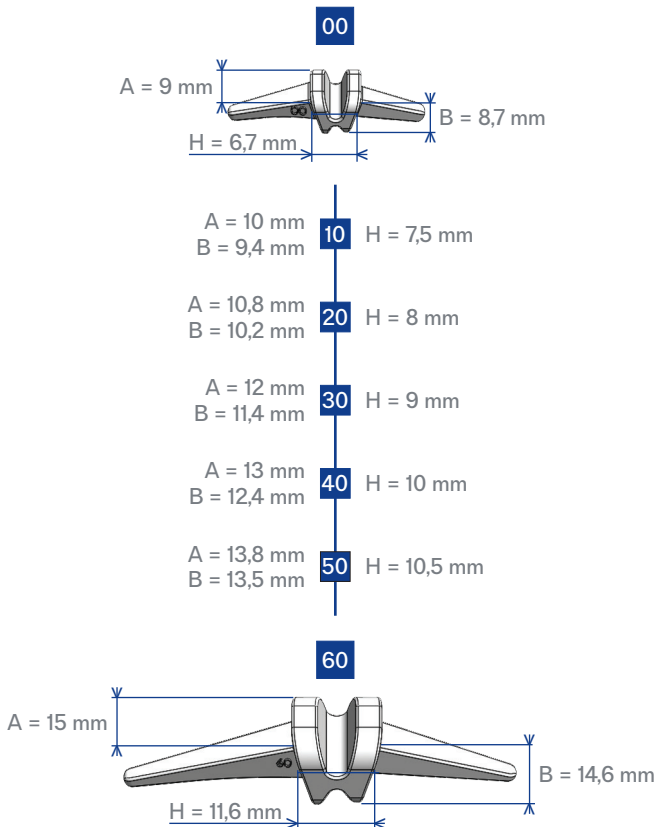




## OPTIMISED & SIMPLIFIED INSTRUMENTATION


### Instruments

2 progressive rasps per indication

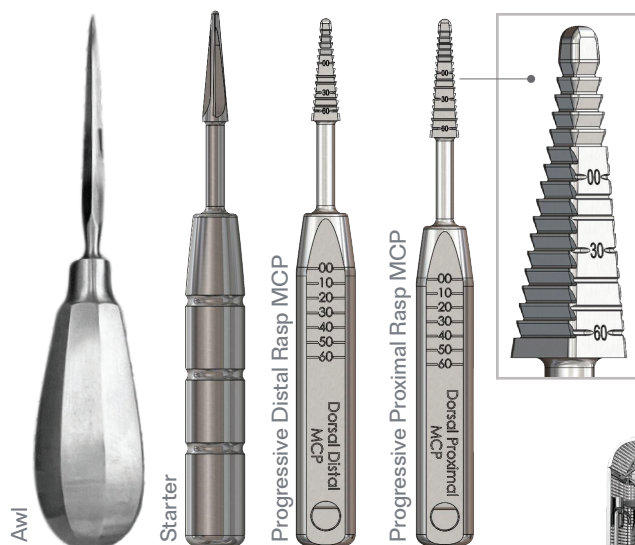


Validated sterilisation cycle 134°/18 min and 134°/3 min according to ANSI/AAMI ST79 and EN ISO 17665


### PIP Instruments

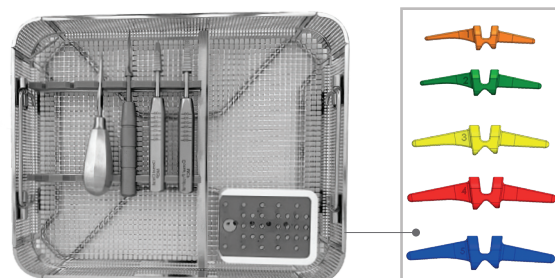

### MCP Instruments

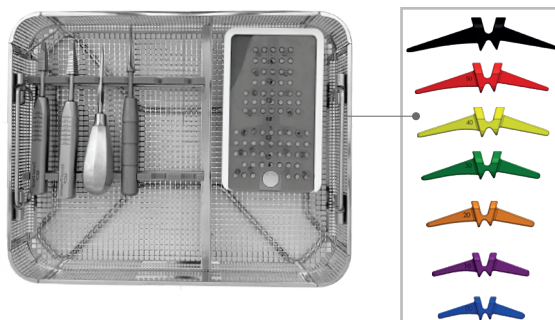

### PIP + MCP Instruments

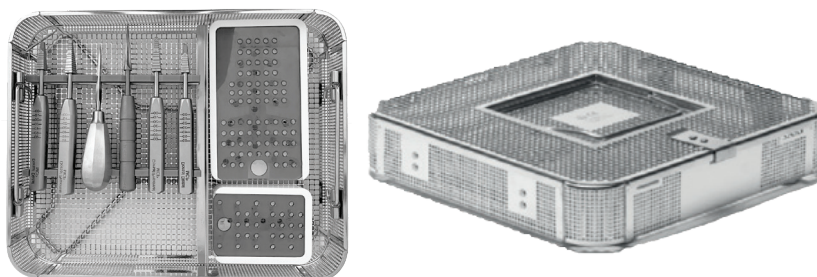

## ORDERING INFORMATION

### PIP Implants

#### Catalog

| Ref        | Designation    | Sizes  |
|------------|----------------|--------|
| 210-P20001 | KeriFlex PIP 1 | Orange |
| 210-P20002 | KeriFlex PIP 2 | Green  |
| 210-P20003 | KeriFlex PIP 3 | Yellow |
| 210-P20004 | KeriFlex PIP 4 | Red    |
| 210-P20005 | KeriFlex PIP 5 | Blue   |

### MCP Implants

#### Catalog

| Ref        | Designation     | Sizes  |
|------------|-----------------|--------|
| 210-P10000 | KeriFlex MCP 00 | Blue   |
| 210-P10010 | KeriFlex MCP 10 | Purple |
| 210-P10020 | KeriFlex MCP 20 | Orange |
| 210-P10030 | KeriFlex MCP 30 | Green  |
| 210-P10040 | KeriFlex MCP 40 | Yellow |
| 210-P10050 | KeriFlex MCP 50 | Red    |
| 210-P10060 | KeriFlex MCP 60 | Black  |



**KeriMedical SA**  
Switzerland (HQ)  
Route des Acacias, 45a  
1227 Geneva - Switzerland  
Tél: +41 58 255 01 30  
Fax: +41 58 255 01 40  
customerservice@kerimedical.com

**KeriMedical SAS**  
France  
ActiPro - 219 rue Laszlo Biro  
74160 Archamps - France  
Tél: +33 (0)4 50 85 36 10  
Fax: +33 (0)4 80 16 05 10  
adv@kerimedical.com

**KeriMedical GmbH**  
Germany/Austria  
Friedrich-Koenig-Str. 3  
55129 Mainz - Germany  
Tel: +49 (0) 6131 617 69 26  
Fax: +49 (0) 6131 617 69 28  
officegermany@kerimedical.com

