



KeriFlex[®]

Silicone Finger Joint Implants

PIP - MCP



Fully committed to **Hand & Wrist** only
www.kerimedical.com

KeriFlex®

Silicone Finger Joint Implants PIP - MCP

Replaces joints affected by rheumatoid arthritis, osteoarthritis or traumatic arthritis


+ COMFORT

+ STABILITY

KeriFlex® implants provide strength, flexibility and durability

**MEDICAL SILICONE
and BIOCOMPATIBLE**

- Anatomic center of rotation is respected
- Hinge design adapted to the joint, limiting conflicts with the extensor mechanism
- Resistance to natural compression stress of the joint

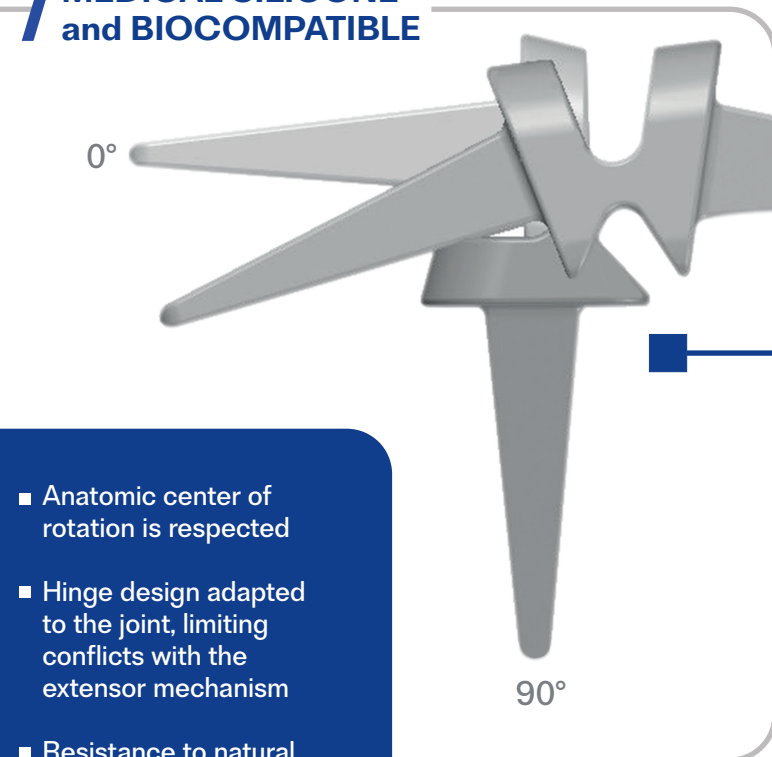


0°

90°



Test lab

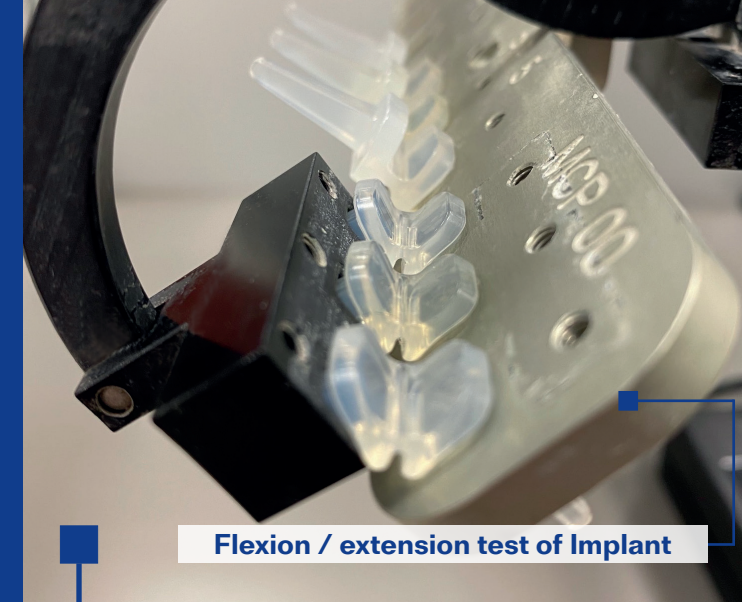

Flexion / extension test of Implant

DESIGN & RESISTANCE VERIFIED

- Full flexion / extension (0° - 90°)
- Verified resistance up to 10 million cycles (test lab)
- Respect of joint angulation at rest : PIP 15° and MCP 30°

The new generation
of Finger Joint Implants

Medio-lateral stability
and joint mobility


High performance
medical silicone
for high mechanical
resistance
to physiological efforts

Complete size range to
meet all anatomical needs

12 SIZES AVAILABLE


PIP

5 sizes of implants

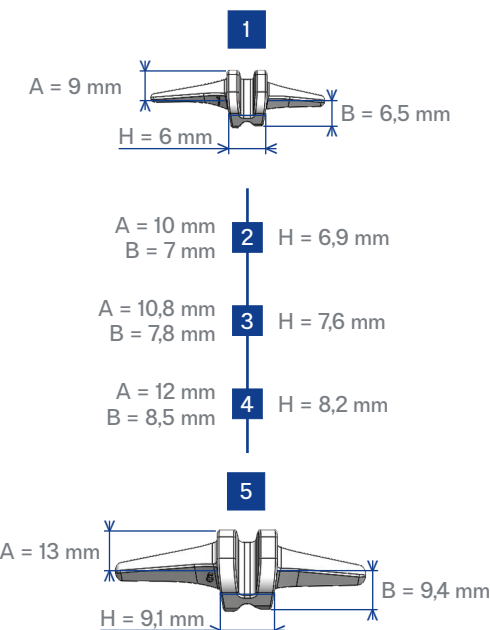

MCP


7 sizes of implants

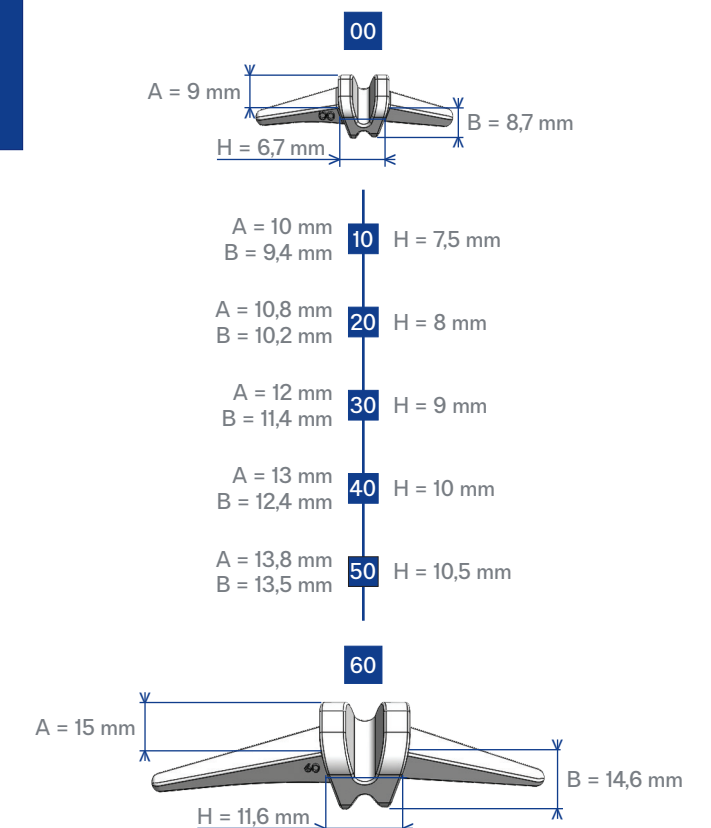

OPTIMISED & SIMPLIFIED INSTRUMENTATION


Instruments

2 progressive rasps per indication

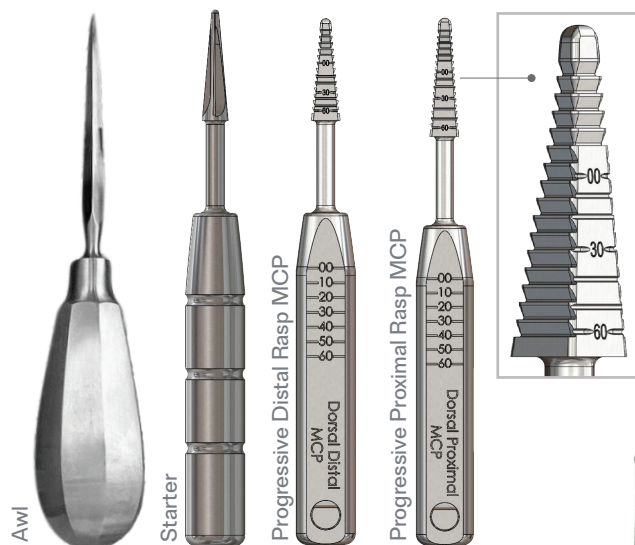


Validated sterilisation cycle 134°/18 min and 134°/3 min according to ANSI/AAMI ST79 and EN ISO 17665

Set PIP

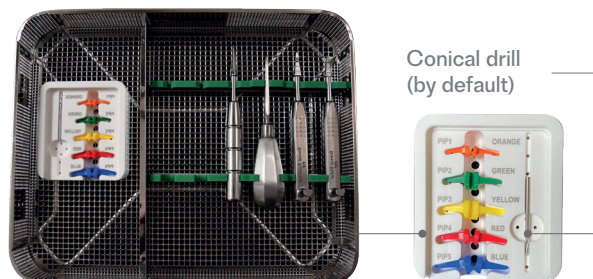


Conical drill (by default)

Set MCP

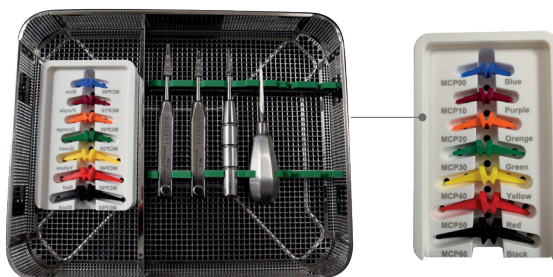

Set PIP + MCP

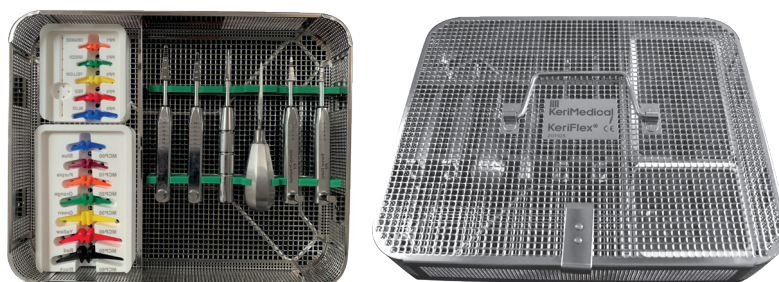

ORDERING INFORMATION

PIP Implants

Catalog

Reference	Designation	Pattern
210-P20001	PIP 1 Finger Joint Implant	Orange
210-P20002	PIP 2 Finger Joint Implant	Green
210-P20003	PIP 3 Finger Joint Implant	Yellow
210-P20004	PIP 4 Finger Joint Implant	Red
210-P20005	PIP 5 Finger Joint Implant	Blue

MCP Implants

Catalog

Reference	Designation	Pattern
210-P10000	MCP 00 Finger Joint Implant	Blue
210-P10010	MCP 10 Finger Joint Implant	Purple
210-P10020	MCP 20 Finger Joint Implant	Orange
210-P10030	MCP 30 Finger Joint Implant	Green
210-P10040	MCP 40 Finger Joint Implant	Yellow
210-P10050	MCP 50 Finger Joint Implant	Red
210-P10060	MCP 60 Finger Joint Implant	Black



KeriMedical SA
Switzerland (Headquarters)
 Route des Acacias, 45a
 1227 Geneva - Switzerland
 Tél: +41 58 255 01 30
 Fax: +41 58 255 01 40
 customerservice@kerimedical.com

KeriMedical SAS
France
 ActiPro - 219 rue Laszlo Biro
 74160 Archamps - France
 Tél: +33(0)4 50 85 36 10
 Fax: +33 (0)4 80 16 05 10
 adv@kerimedical.com

Distributed by:

